

A-F-T
Automatic Flexion




Hill Laboratories
COMPANY

AIR AFT

Auto-Flex with Air-Drops

Crisp, smooth and impressively fast.

Equipped with the Hill Air-Drop System*, the AirAFT allows doctors to select one or up to four drops simultaneously from the convenient, side-mounted, control panel. Cushions are activated instantly by tapping one of the foot strips.

Shown in Cinder Upholstery with Black Base.



** Additional options shown including Raised Headpiece, Gripper Bar variation, Foot Strips and Timer.*

Select one or up to four drops simultaneously from the Control Panel. Optional dual-sided control box available

** The Hill Air-Drop System, with air-powered features, is only available on the AirAFT.*

STANDARD FEATURES

- Variable speed motorized flexion
- Air activated cervical, thoracic, lumbar and pelvic drops
- Air cushioned thoracic breakaway and thoracic release
- Tilting Headpiece, 30° negative and positive tilt
- Pelvic section has lateral flexion and can be locked at any angle up to 21°
- Removable tiller bar at foot-end for lateral movement
- Slide-out ankle support with ankle harnesses
- Adjustable patient gripper bar
- Patient Shut-Off Switch
- Electric height controlled by Dual Rocker Foot Pedals
- Width 24" / Height 21.5" to 29", higher heights available
- Foam Top - 2 1/2" Medium (standard), Soft or Firm
- Choice of 20 Colors for the top cushions
- Beige, Gray or Black Base, enclosed for safety
- One-year warranty on all parts and labor

HILL AFT™

Motorized Automatic Flexion

The Hill Automatic Flexion Table is perfect for doctors who want full-featured automatic flexion without the need for manual flexion. The AFT is surprisingly affordable for a table that features variable speed flexion, electric adjustable height and an impressive list of other standard features. Make the AFT even more versatile by adding up to four manual drops and several other options to customize the AFT for your specific needs.



Pictured with several options

Shown in Navy Upholstery with Beige Base.

▼ STANDARD FEATURES

- Variable speed motorized flexion
- Tilting Headpiece, 30° negative and positive tilt
- Thoracic release and lock
- Pelvic section has lateral flexion and can be locked at any angle up to 21°
- Removable tiller bar at foot-end for lateral movement
- Slide-out ankle support with ankle harnesses
- Adjustable patient gripper bar
- Patient Shut-Off Switch
- Electric height controlled by Dual Rocker Foot Pedals
- Width 24" / Height 21.5" to 29", higher heights available
- Foam Top - 2 1/2" Medium (standard), Soft or Firm
- Choice of 20 Colors for the top cushions
- Beige, Gray or Black Base, enclosed for safety
- One-year warranty on all parts and labor



The AFT has adjustable height 21.5" to 29"



Slide-Out Ankle Support

AFT STANDARDS

All of the following features are standard on the AFT

AUTOMATIC MOTORIZED FLEXION

Variable Speed Flexion and much more.

The Hill AFT has a Flexion Control Panel that is conveniently located at the foot-end of the table which allows you to set the flexion speed from 1 to 28 rpm (an optional 30 minute timer can also be added). A selector switch allows you to stop the table in the horizontal or flexed position. In the flexed position, with the table at the lowest height, patients can easily enter or exit the table from the foot end.

The pelvic section of the AFT moves laterally and can be easily maneuvered during flexion using the Tiller Bar. It can also be locked laterally to either side.

The Hill AFT is shown here with optional manual drops including Dual-Drop Cervical, Thoracic, Lumbar and Pelvic.



Shown in Cinder and Black



TILTING HEADPIECE

Hand lever releases the headpiece to lock at any angle from 30° positive to 30° negative.



TILLER BAR

Inserts for using manual lateral movement or circumduction during a flexion treatment.



THORACIC RELEASE

Thoracic Release locks at any angle down to 4 1/4" and up (above horizontal) to 1 1/2".



PATIENT GRIPPER BAR

Adjusts in length to suit patient size and comfort during flexion.



ROCKER FOOT PEDAL

Mounted on both sides of the table base for easy height adjustment.



PATIENT SHUT-OFF

Equipped with a patient shut-off switch located near the left armrest.

AFT OPTIONS

Made-to-Order

▶ MANUAL OR AIR DROPS


The AFT can be equipped with up to 4 manual drops including cervical (straight-down and forward-motion), thoracic, lumbar and pelvic. The AirAFT comes standard with all 4 air drops.

▶ THORACIC SPRING BREAKAWAY

For a thoracic recoil movement, a fixed-tension gas shock can be added to the AFT. This feature can be combined with a thoracic or lumbar drop. The AirAFT comes standard with an air-controlled breakaway and pressure pedal.


▶ HEADPIECES

Add functionality to your headpiece by combining several headpiece options including Dual-Drop, Raised, Flexion, Axial-Motion, Lateral and Translation (CBP). The Manual AFT comes standard with a Tilt headpiece. The AirAFT comes with a Tilt and Dual-Drop headpiece.

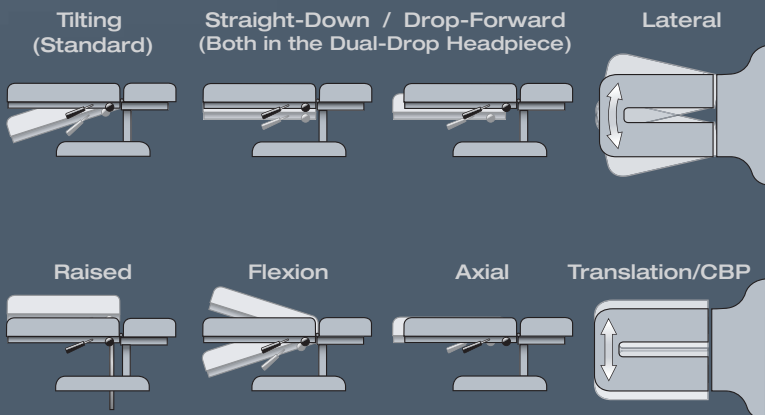

▶ MANUAL CERVICAL FLEXION

The Flex-Headpiece (shown right) comes with the raised headpiece, occipital strap and removable tiller bar. The headpiece unlocks and a gas-spring gently compresses and recoils as you control the flexion, lateral and long-axis movement.

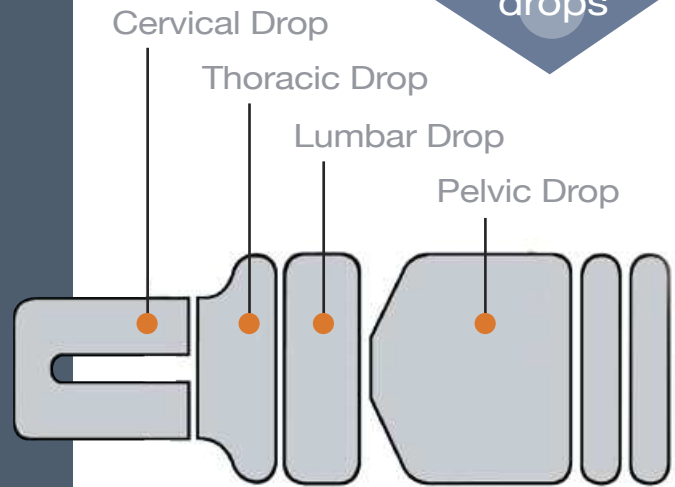

Thoracic/Abdominal Breakaway



DUAL FOOT STRIPS

With a tap of the foot, you can activate the drops from either side of the table.



CASTERS

Mount on each end of the AFT base to permit easy movement of the table.



SPECIFICATIONS

Electric actuator 60 Hz, 110 volts, 4.8 amps./AirAFT - 10 amps. (220-50 cycle also available)

Lifting capacity - 450 lbs.

Length - 5' 5" (extends to 6' 2"), with cervical drop - 5' 7" to 6' 4" • Width - 24"

Range of height; 21 1/2" to 29"; higher heights available

Medium density foam top; firm and soft density upon request

Choice of colors for upholstered top (see below)

Beige, gray or black telescoping base skirting constructed of fire-rated thermo-plastic

Powder-coated internal and external metal parts

Shipping weight 250-300 lbs.

One-year warranty on all parts and labor

VERY AFFORDABLE


Finance your Hill Laboratories table through our convenient and affordable lease/purchase plan.

Visit HillLabs.com to calculate leasing options or call to speak personally with one of our representatives who will be happy to discuss payment options.

AVAILABLE COLORS



Below - Finesse® Colors. Upgraded material; foam-backed and puncture-resistant.

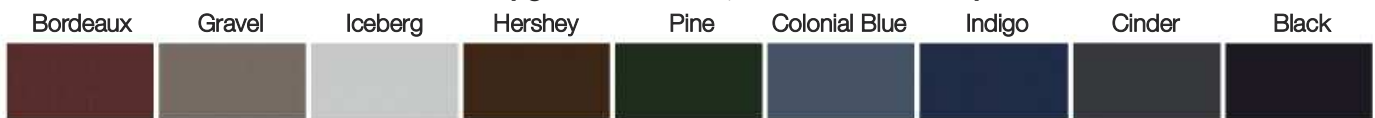

Base Colors



Upholster your table in durable Boltaflex® with Prefixx™, scientifically formulated to give superior protection from staining and abrasions. Finesse and custom colors are available at additional cost. Shades may vary from actual material.

Color choices may change without notice.



P.O. Box 2028, 3 N. Bacton Hill Road, Frazer, PA 19355 • Toll-Free 1-877-445-5020
610-644-2867 • fax 610-647-6297 • www.HillLabs.com