
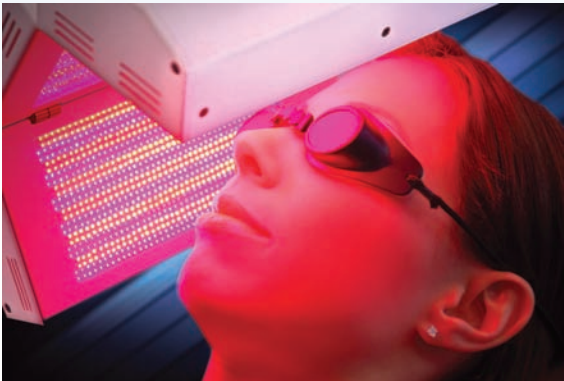


AKLARUS

DUAL WAVELENGTH
PHOTOTHERAPY SYSTEM



The **AKLARUS Blue and Red** Combination System is an effective, drug-free alternative for treating acne and photodamaged skin. The non-invasive Aklarus treatment has been shown to significantly reduce inflammatory acne lesions. The narrow band blue light destroys the P. Acne bacteria while the Red light penetrates deeper stimulating cellular mechanisms responsible for tissue repair and regeneration.



Aklarus Anti-Aging **Infrared and Red** Combination System is used to treat wrinkles and photodamaged skin by stimulating collagen synthesis and increasing the amount of dermal elastic fibers.

Aklarus Wavelength Specifications: Blue, Red and Infrared

Blue 420nm
10W+/- 3W Total Output
Standard dose 26J/ cm²
in 20 Minutes

Red 628nm
20W +/- 3W
Standard dose 52J/cm²
in 20 Minutes

Infrared 880nm
16W +/- 3W
Standard dose 42 J/cm²
in 20 Minutes



BEFORE & AFTER

AFTER 4 WEEKS / 8 TREATMENTS

These results have not been retouched in any way.

