

THE HILL ADJUSTABLE MODEL HA90V

featuring electric tilt, electric adjustable height and manual lift-back



For:

- Endovenous Treatment
- Vascular Diagnostic Testing
- Echocardiogram Testing

Basic Price
\$3495

The vascular exam table Model HA90V is a versatile table offering a sturdy welded-steel base construction and smooth electric elevation.

Standard Features:


- Electric Adjustable Height Base (range 24" to 34" or 26" to 39") controlled by foot pedal
- Electric Tilting Table Top - (head and foot end tilt 20°) also controlled by foot pedal
- Manual Lift Back with switch-activated gas-spring
- Top thickness - 2 1/2" medium density foam . Top size - 27" x 5'11"
- Lifting Capacity - 500 lbs. • Approximate weight - 225 lbs.
- Choice of 20 upholstery colors with gray or beige base

Optional Features:

- | | | |
|-----|---|---------------|
| 98 | 30" Extra Width |\$50.00 |
| 88 | 32" Width |\$100.00 |
| 93 | Caster System for mobility. Engaged with foot lever |\$150.00 |
| 60 | Fold-down Echo Window (* see note below) |\$75.00 |
| 94V | Electric Lift-Back |\$495.00 |
| 48 | Seat Belt |\$40.00 |
| 49 | Foot Board |\$150.00 |
| 50 | Protractor |\$40.00 |
| 140 | Upgrade Upholstery, Puncture-Resistant Tricot Finesse | ...\$225.00 |
| 96 | Paper Roll Bracket |\$25.00 |

* The echo-window features an easy-access 8"x11" fold-down section designed for left lateral position echo cardiology (only available with 32" width)

O P T I O N S

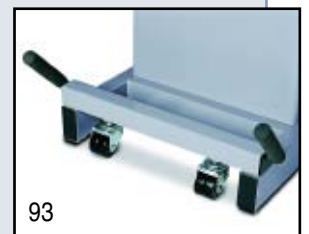

Hill Laboratories
COMPANY