

# The Hill Adjustable

## Air-Drop HA90C

*The first AFFORDABLE Chiropractic Table  
with pneumatic drops!*



**I**f you have used a manual drop table for any length of time, then you'll appreciate the impressive speed and ease of treating with the Air-Drop HA90C. Equipped with the Hill Air-Drop System, the HA90C allows doctors to select one or up to four drops simultaneously from the convenient, side-mounted, control panel. Once a section has been dropped, it is instantly reset with the simple touch of a foot pedal.


All of the pneumatic components on the Air-Drop HA90C are supplied by a large air-storage tank and mini-compressor which are neatly enclosed within the table's base skirting. The compressor runs quietly and only periodically (approximately 1 minute) to replenish the air-tank.

The table features an Air-Breakaway controlled by a foot pedal. The pedal increases or decreases the air-spring pressure in the thoracic and lumbar sec-


tions providing a controlled recoil action. It can also be locked in any position down to 5½" for treating patients with a large abdominal region.

The drops included on the HA90C are Cervical and Pelvic, plus a straight-down Thoracic and Lumbar drop. The Cervical Drop can be equipped with a Forward-Motion or a Straight-Down action. An optional Raised Headpiece can also be added.

Like all Hill Adjustables, the Air-Drop HA90C has electric variable-height. Height is adjusted with a rocker foot pedal that is mounted to the base and can be accessed from either side of the table. An operating range from 23½" to 32" allows the doctor to work at a comfortable treatment level. Tables are available in a choice of twenty upholstery colors with a base of either beige or gray. All Hill Tables come with a one-year warranty on all parts and labor.

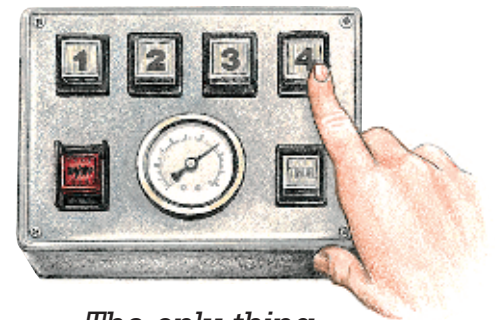


*The only thing  
easier than  
selecting the drops...*

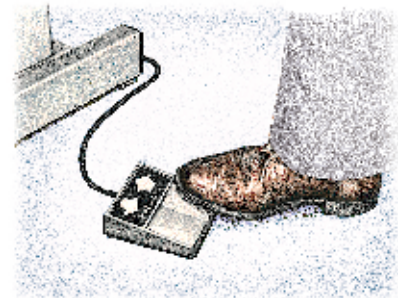

*is cocking them!*

*- Please see other side for more important information -*


# Standard Features

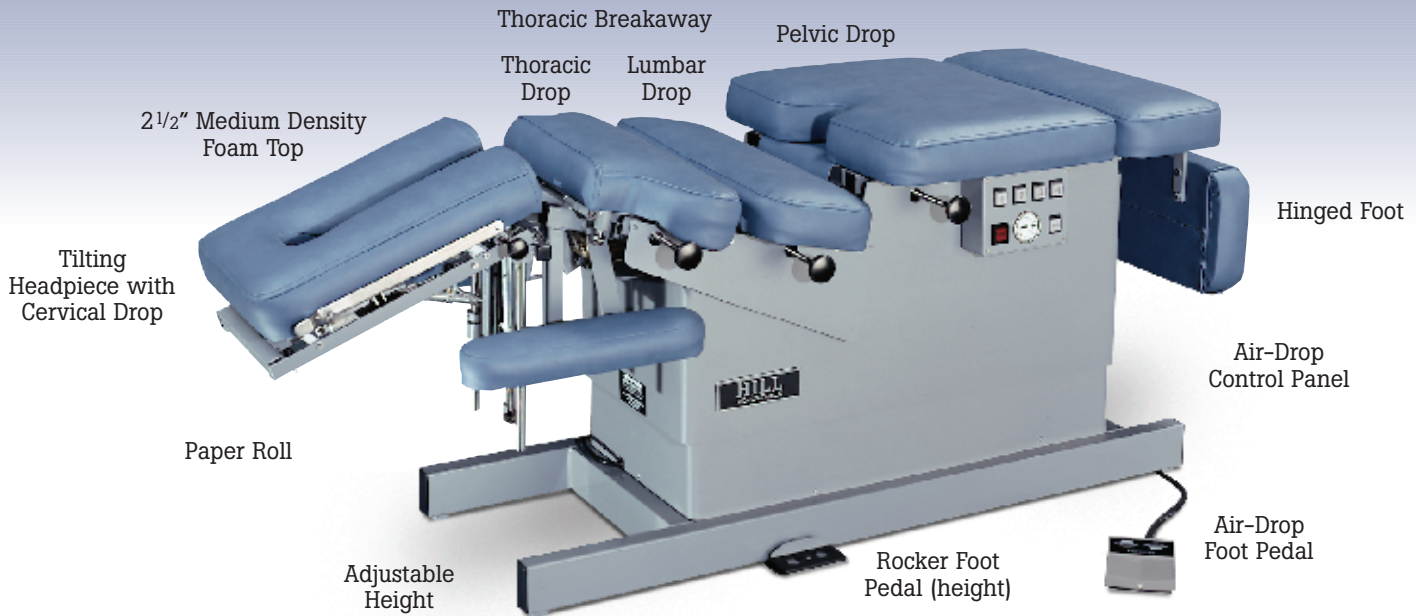


# Optional Features




Tilting Headpiece with forward-motion or straight-down drop

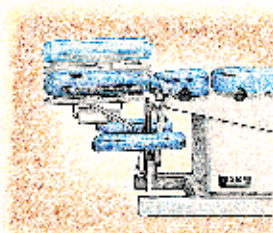

Raised Tilting Headpiece with forward-motion or straight-down drop

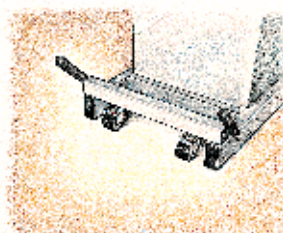

Caster System

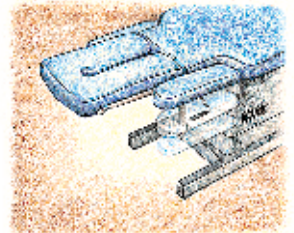

Adjustable Arm Rests

# Specifications

- Powered by a thermally protected electric actuator - 60 cycle, 110 volts (220-50 cycle additional charge) and a thermally protected low-noise air-compressor. Total amps - 7.0 • Range of height on HA90C, from 23 1/2" to 32"
  - Lifting capacity - 450 lbs. • Standard width, 24"; optional widths, 22", 27", or 30"
- Table top length - 6'1"; with hinged foot down - 5'4" • 2 1/2" medium density foam
  - Beige or gray telescoping base constructed of fire-rated acrylic thermo-plastic
    - Choice of 20 colors for upholstered top. Special colors upon request
    - Table shipping weight approximately 250 lbs.
    - One-year warranty on all parts and labor.

Air-Drop convenience, great Hill styling, sturdy craftsmanship, and surprising affordability make the Air-Drop HA90C another great product in the Hill tradition of value.

*Hill Laboratories Company is a family-owned business serving the chiropractic profession for over 55 years.*

  
**Hill Laboratories**  
**COMPANY**

P.O. Box 2028, 3 Bacton Hill Road, Frazer, PA 19355 • 610-644-2867 • fax 610-647-6297