

# ElectroLight™

## Infrared Light Therapy Device



**Only \$1595**

*With 4 Infrared Applicators and  
5 Packs of Dual-Sided Electrodes.*

### Add Infrared light to your stimulation unit for combination therapy.



The Electro-Light is an infrared therapy device that can be used in conjunction with a separate muscle stimulation machine. The Electro-Light uses a patented transparent electrode called the Dual-Sided Electrode™. The Dual-Sided Electrode is clear so infrared light can pass through the electrode at the same time as muscle stimulation. The Dual-Sided Electrodes can be placed on a patient similar to a standard electrode. The Dual-Sided Electrodes are sticky on both sides so the Electro-Light applicators can stick to the electrode or a strap can secure the applicator to the top-side of the electrode. Any stimulation machine can be used in combination with the Electro-Light.

#### Electro-Light Features:

- Electro-Light Console with Adjustable Time and Frequency Settings
- 4 Infrared (880nm) and Red(625nm) Thermal Applicators with 340mW of Power Each
- 5 Packs of Dual Sided Electrodes (2" x 2" square)
- Patient Cut Off Switch
- FDA Clearance

#### Additional Therapy Add-Ons:

- Rolling Cart - \$100
- Deluxe Rolling Cart with Metal Basket - \$130
- Dual Sided Electrodes 2" x 2" Square (4)(Quantity discounts Available) - \$3.50
- 48" x 2" Velcro Strap - \$9.00

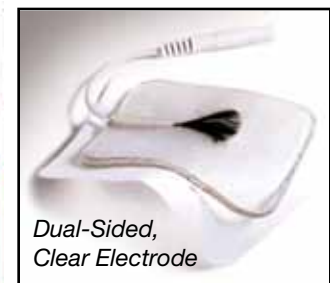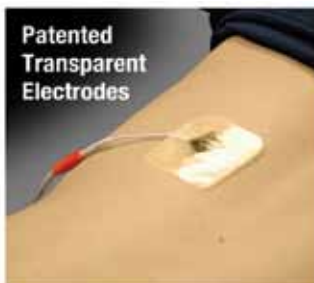

**1-877-445-5020 • [www.HillTherapy.com](http://www.HillTherapy.com)**  
3 N. Bacton Hill Road, Frazer, PA 19355  
[info@HillTherapy.com](mailto:info@HillTherapy.com)

 **Hill Therapeutics**  
A Division of Hill Laboratories